

Argon Purity

Materials Qualification for long drift Liquid Argon TPC

Effect of water on electron drift lifetime

**Plans Material Tests, purging and H₂O studies,
Purity in a large tank without evacuation (LAPD) – Brian Rebel
MicroBooNE**


**Acknowledge enormous and continuing contribution of ICARUS and
other European programs**

Argon Source, Materials Test System, and Electronics Test System



Develop purification, check detector materials, TPC for electronics development

Schematic of Materials Test Stand

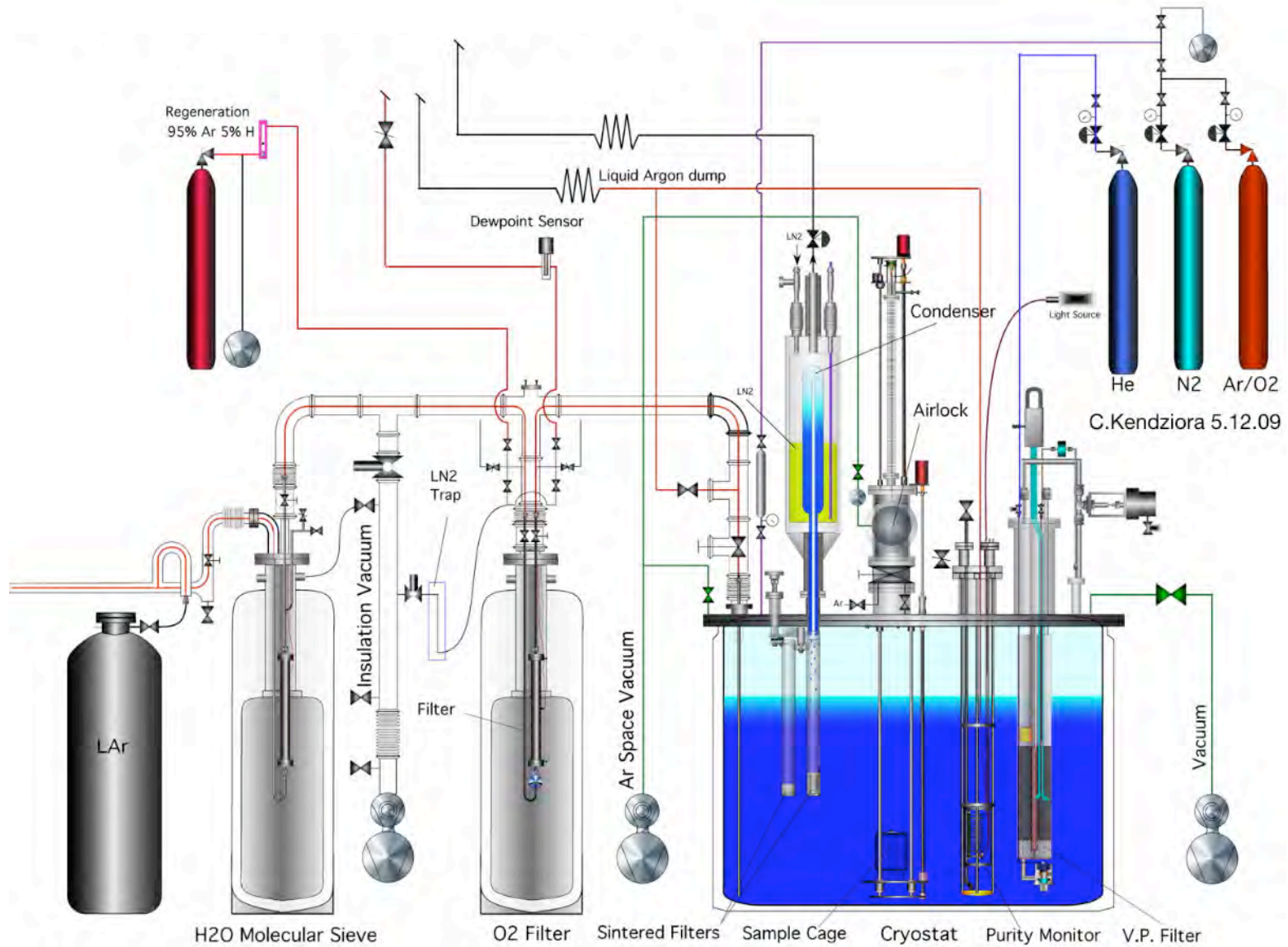

Materials Test Stand





Features of Materials Test Stand

- Can insert materials into known clean argon
- Can insert materials after purging only or after pumping on them.
- Can position materials into liquid and into ullage with range of temperatures.
- Can insert known amounts of contaminant gases
- Nitrogen-based condenser can maintain liquid for long (weeks) studies
- Internal filter-pump can remove contamination introduced by materials – 2hr cycle
- Sample points at Argon Source, after single-pass filters, in cryostat gas and liquid

Measurement Features of Materials Test Stand

- Measure electron drift lifetime with ICARUS style purity monitor
- Measure Oxygen (0.3 ppb sensitivity) with oxygen meter (Delta-F & Tiger Optics)
- Measure H₂O in gas (0.3 ppb sensitivity) with water meter (Tiger Optics)
- Cryogenic data, Lifetime Data, analytic instrumentation data in single data-base


Insertion of a material sample into the airlock (argonlock) basket

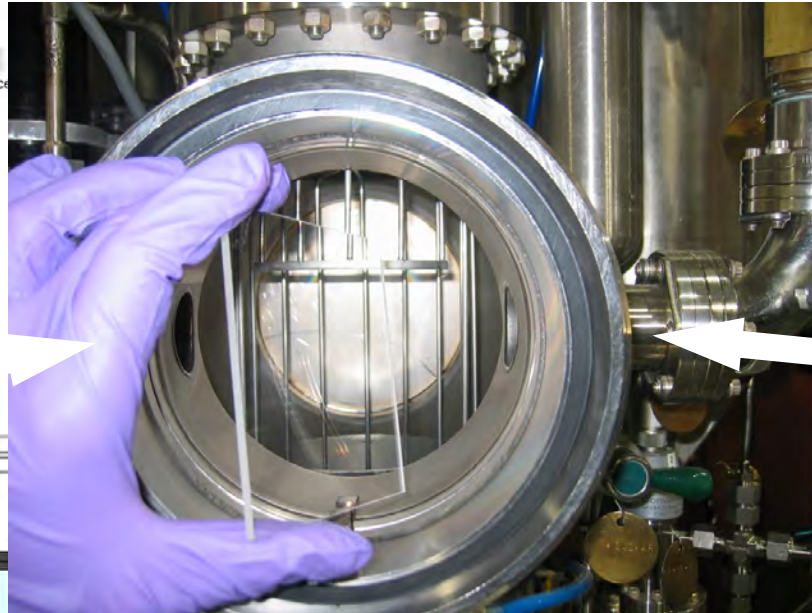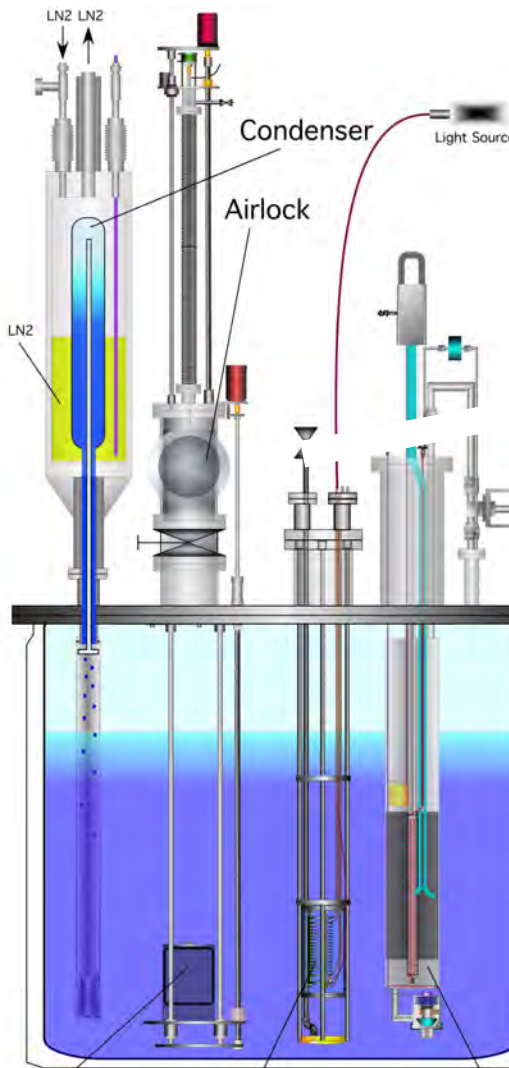



Sample volume 10 cm x 10 cm x 10 cm
Argonlock can be purged with
external argon, cryostat argon, and/or
evacuated

Sample Cage Purity Monitor Scrubber Filter


Cryogenics Control Screen

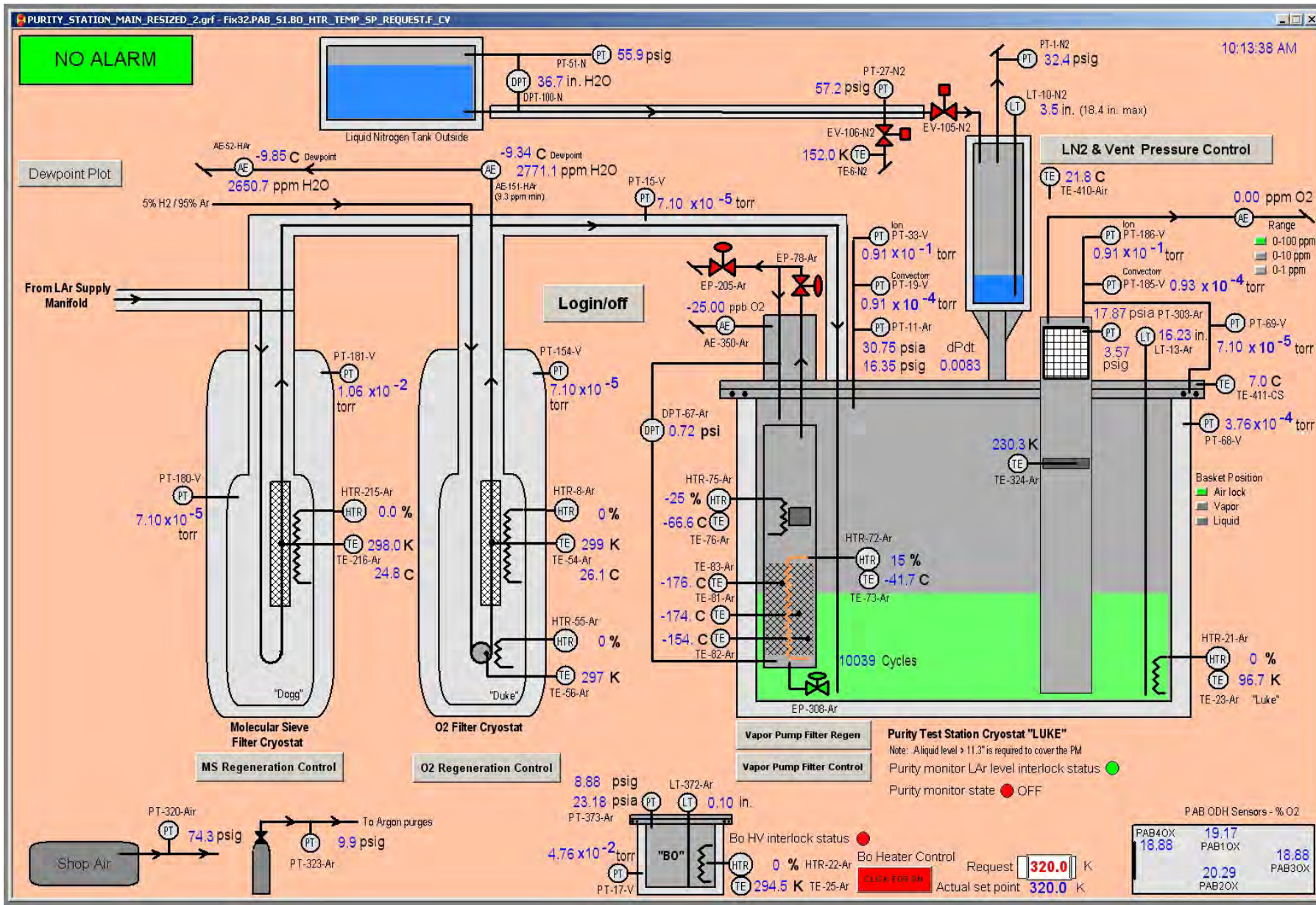

Electron Drift Lifetime Monitor

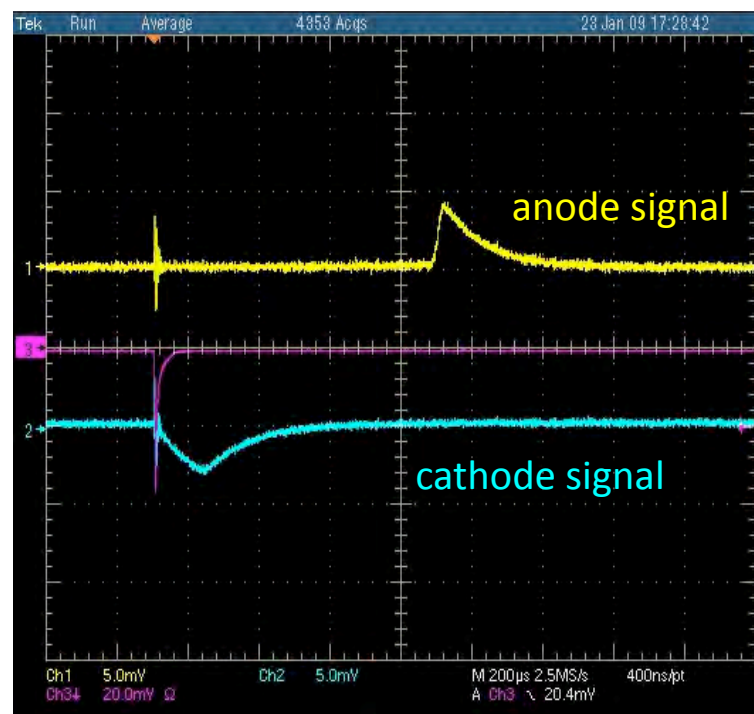
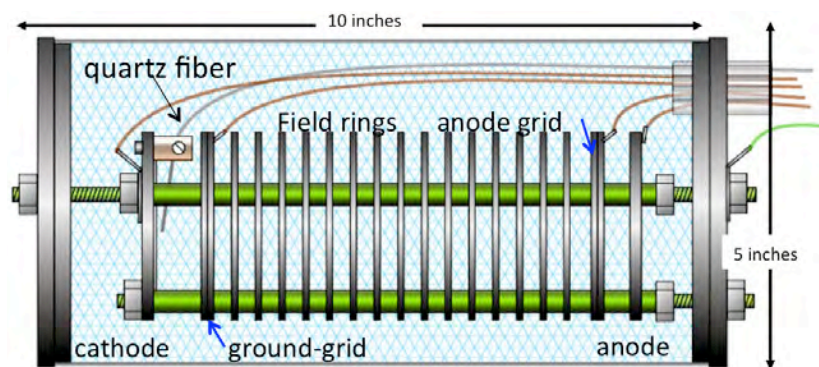


PrM scope signal

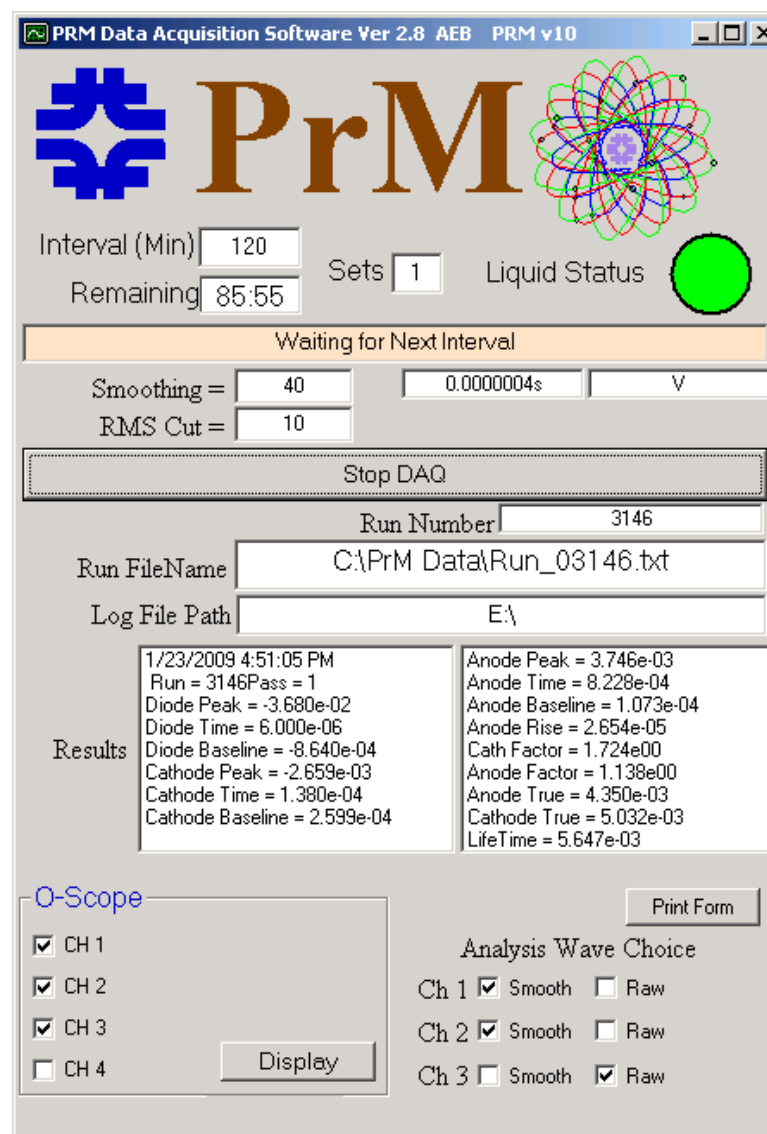

PrM automation software

An example set of measurements with the Materials Test System

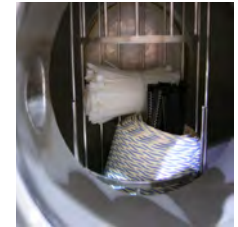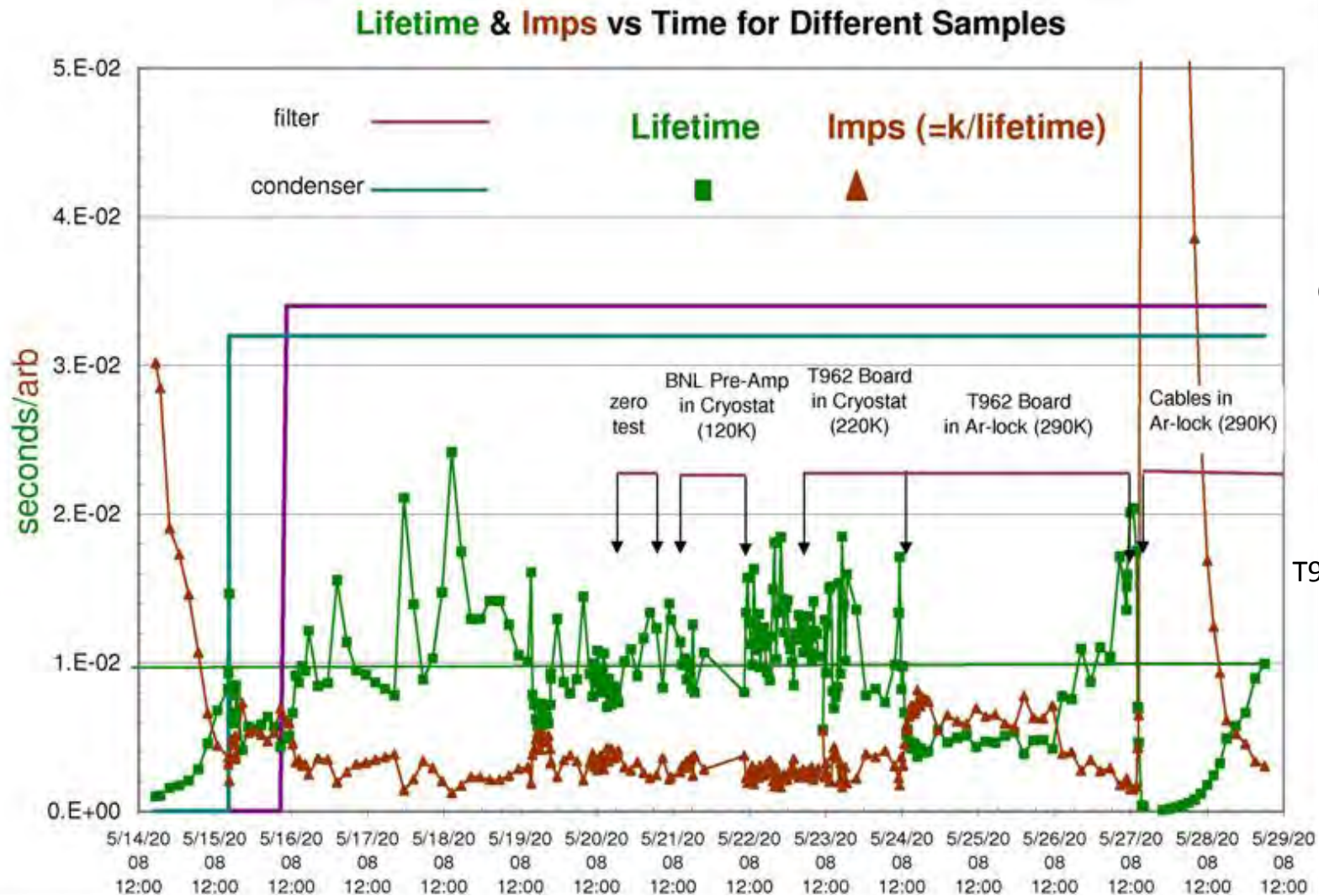

Cables & Cable ties



T962 DecouplingBoard



Cold pre-amp

Learning how to do what has been done by others

(cryogenics, purification, purity monitoring - all are now designed & built in the US)

New stuff - our own filter systems,
material test systems,
the effect of H₂O,

FERMILAB-TM-2384-E: efficiency of slow purging to remove
atmosphere to ppm levels

A regenerable filter for liquid argon purification

A. Curioni ^b, B.T. Fleming ^b, W. Jaskierny ^a, C. Kendziora ^a, J. Krider ^a, S. Pordes ^a, M. Soderberg ^b,
J. Spitz ^{b,*}, T. Tope ^a, T. Wongjirad ^b

NIM-A 605:306-311,2009.

^a Particle Physics Division, Fermi National Accelerator Laboratory, Chicago, IL, USA

^b Department of Physics, Yale University, New Haven, CT, USA

A system to test the effect of materials on electron drift lifetime in liquid argon and the effect of water

R. Andrews, W. Jaskierny, H. Jöstlein, C. Kendziora, S. Pordes ^{*}, T. Tope NIM-A 608:251-258,2009.

Particle Physics Division, Fermi National Accelerator Laboratory, Batavia, IL 60510, USA

Some Materials Tested

Material	Sample Surface Area (cm ²)	Effect of Material on Electron Drift Lifetime (LT)			Comments
		94 K liquid	≈120 K vapor	≈225 K Vapor	
Red-X Corona Dope ^a	100	None	None	LT Reduced from 8 to 1 ms; recovery observed.	H ₂ O concentration not monitored.
Deactivated Rosin Flux ^b	200	None	Not Tested	LT reduced from 8 to 1.5 ms recovery observed	H ₂ O concentration not monitored.
FR4	1000	None	Not Tested	LT reduced from 8 to <1 ms	Outgassed enough H ₂ O at 225 K to saturate sintered metal return.
Taconic ^c	600	None	Not Tested	LT reduced.	Sample outgases water at 225 K.
Hitachi BE 67G ^d	300	None	Not Tested	LT reduced; recovery observed	Sample outgases water at 225K; outgassing reduced over time.
TacPreg ^e	200	None	None	LT reduced; recovery observed	Sample outgases water at 225 K; outgassing reduced over time.
FR4, y-plane wire endpoint for uBooNE	225	None	None	LT reduced from 8 to 3 ms	Sample outgases water at 225 K.
FR4, y-plane wire endpoint for uBooNE	225	None	None	None	Sample was evacuated in airlock prior to testing
FR4, y-plane wire cover for uBooNE	225	None	None	None	Sample was evacuated in airlock prior to testing
Devcon 5-min epoxy	100	None	None	LT reduced from 10 to 6 ms; some recovery observed	Sample outgases water at 225 K.

– no effect on lifetime when material is in Liquid

In summary:


We have a set of tools, hardware and software, to measure effect of materials on electron drift lifetime.

We plan to carry these tools throughout the liquid argon program

Materials Test Stand Program : Test all candidate detector components.

The effect of H2O


FR-4 based circuit board – from Argonlock *with evacuation*

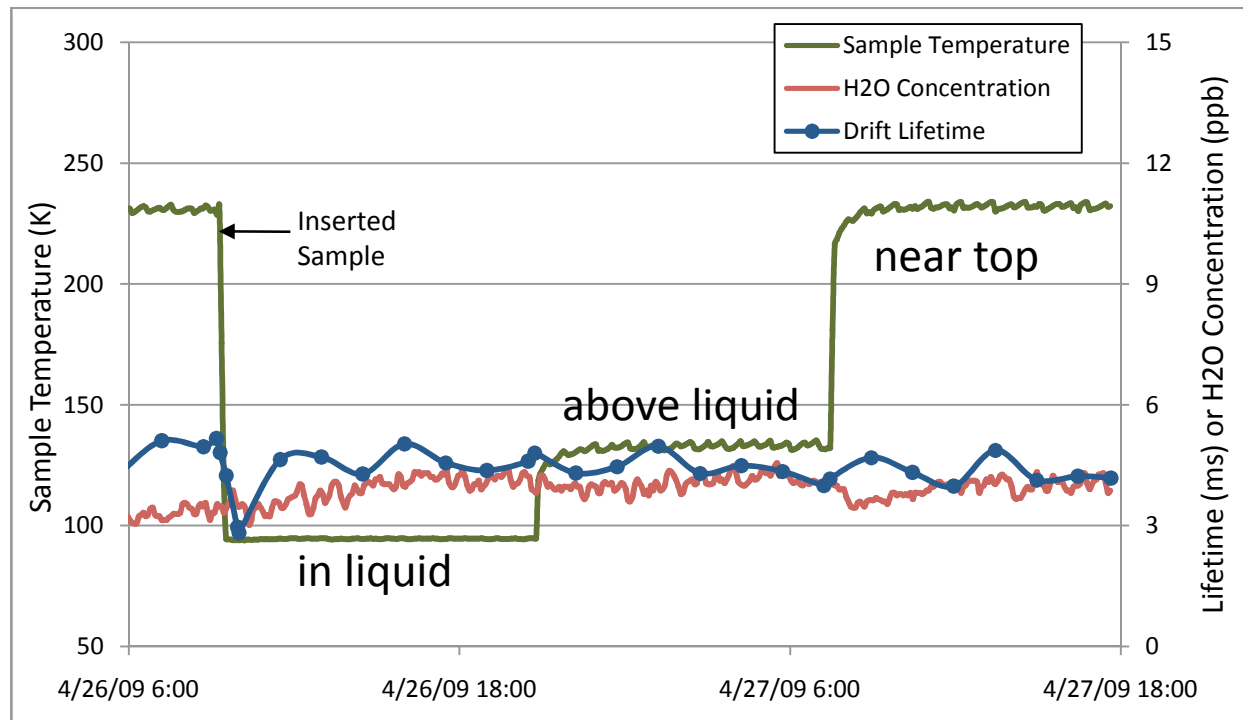

Little change in H2O reading and little change in lifetime

The effect of H2O

FR-4 based circuit board – from Argonlock *with purging only*

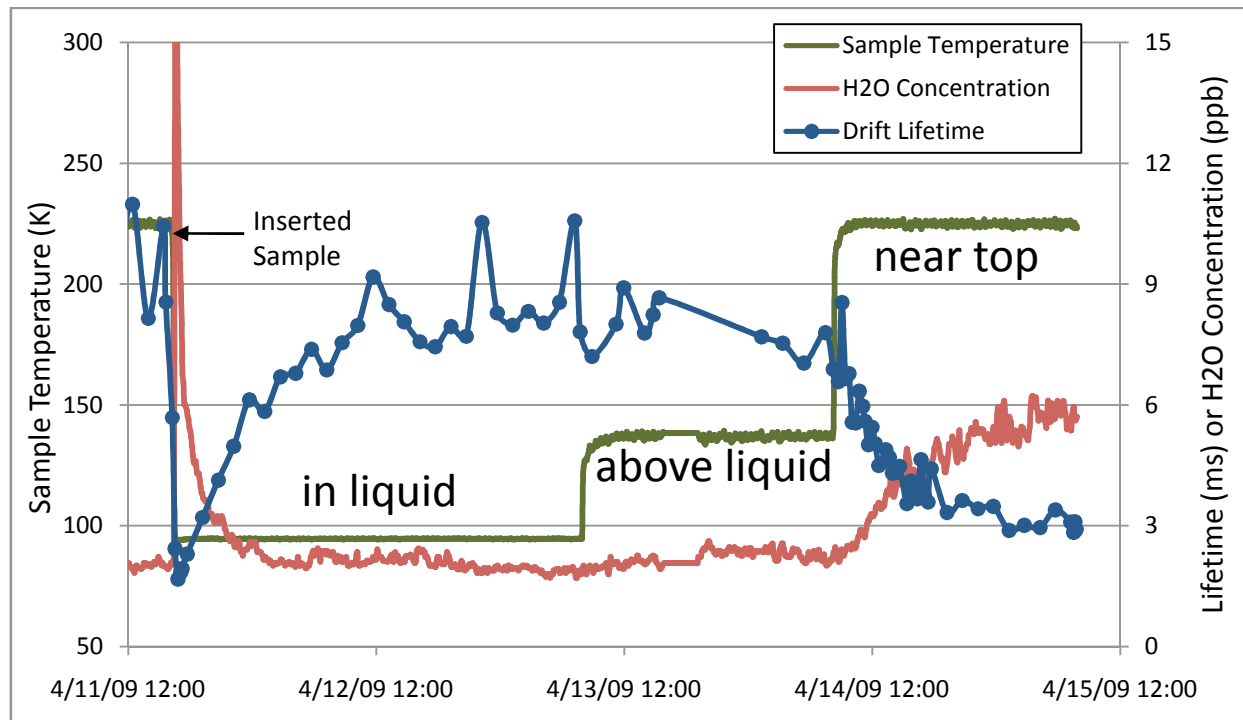

Significant change in H2O reading and significant reduction in lifetime

Also – note initial lifetime better – and H2O reading lower

The effect of H2O

Material	Sample Surface Area (cm ²)	Effect of Material on Electron Drift Lifetime (LT)			Comments
		94 K liquid	≈120 K vapor	≈225 K Vapor	
Red-X Corona Dope ^a	100	None	None	LT Reduced from 8 to 1 ms; recovery observed.	H ₂ O concentration not monitored.
Deactivated Rosin Flux ^b	200	None	Not Tested	LT reduced from 8 to 1.5 ms recovery observed	H ₂ O concentration not monitored.
FR4	1000	None	Not Tested	LT reduced from 8 to <1 ms	Outgassed enough H ₂ O at 225 K to saturate sintered metal return.
Taconic ^c	600	None	Not Tested	LT reduced.	Sample outgases water at 225 K.
Hitachi BE 67G ^d	300	None	Not Tested	LT reduced; recovery observed	Sample outgases water at 225K; outgassing reduced over time.
TacPreg ^e	200	None	None	LT reduced; recovery observed	Sample outgases water at 225 K; outgassing reduced over time.
FR4, y-plane wire endpoint for uBooNE	225	None	None	LT reduced from 8 to 3 ms	Sample outgases water at 225 K.
FR4, y-plane wire endpoint for uBooNE	225	None	None	None	Sample was evacuated in airlock prior to testing
FR4, y-plane wire cover for uBooNE	225	None	None	None	Sample was evacuated in airlock prior to testing
Devcon 5-min epoxy	100	None	None	LT reduced from 10 to 6 ms; some recovery observed	Sample outgases water at 225 K.

Lifetime deteriorates when 'warm' materials out-gas water

The effect of H₂O

In general: Lifetime inversely proportional to H₂O concentration in gas

Constant of proportionality independent of material

– so H₂O is not just a marker and is probably the prime contaminant

Calculations of effective concentration in liquid suggest H₂O as bad O₂

Lessons Learned

Condensate must be filtered before return


Need to understand effectiveness of purging in removing water

Filter Materials

We make and regenerate our own filters using material from BASF (originally Engelhard) for O₂ and Molecular Sieve from Aldrich

We have the capability to use other materials but have not pursued that (yet?)

We are interested in the capacity of filter materials before 'breakthrough' - the capacity seems much less than one would calculate from the quoted area of copper.

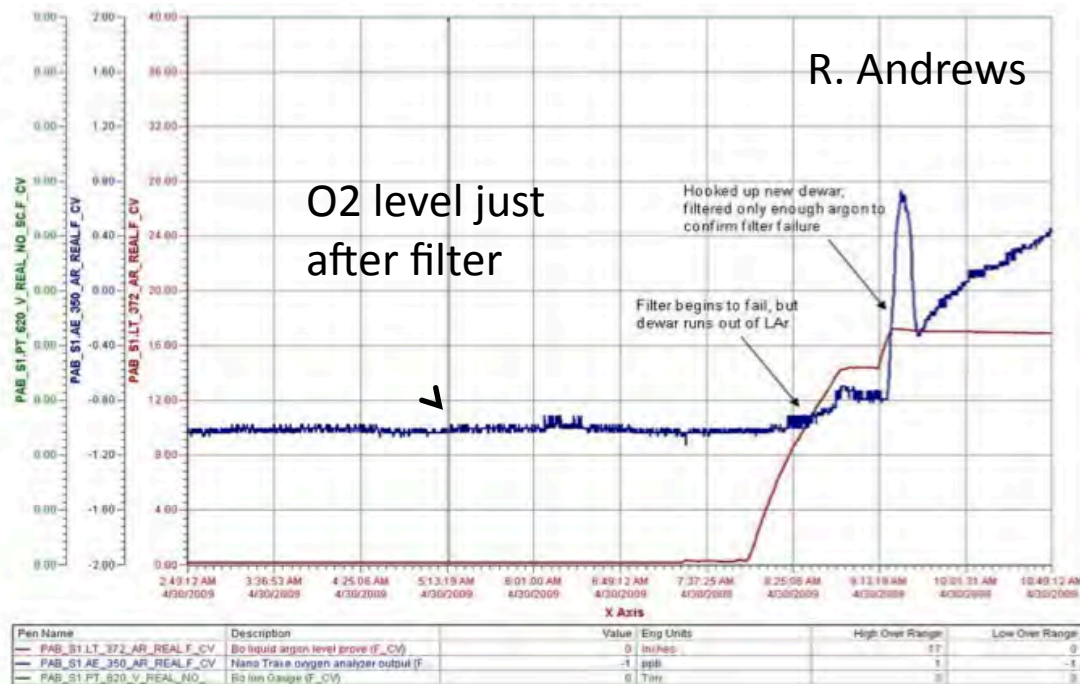

Plot of O₂ level after ~ 0.5 g O₂ through filter. Note the 'breakthrough' style of behavior.

Plans in the Purity Business at PAB

Run materials in the Materials Test Stand

Measure H₂O out-gassing rates of materials

- can infer from the MTS data,
- make a dedicated apparatus (thinking how to)


Measure effectiveness of purging vs time and temperature

Measure Capacity of Filter Materials and parameters thereof
(eg flow rate dependence)

Install Condenser with liquid and gas phase filtration capability
to test effectiveness of each.

**Condenser, filter and
vapor pump all in one.**


(under assembly at
present)

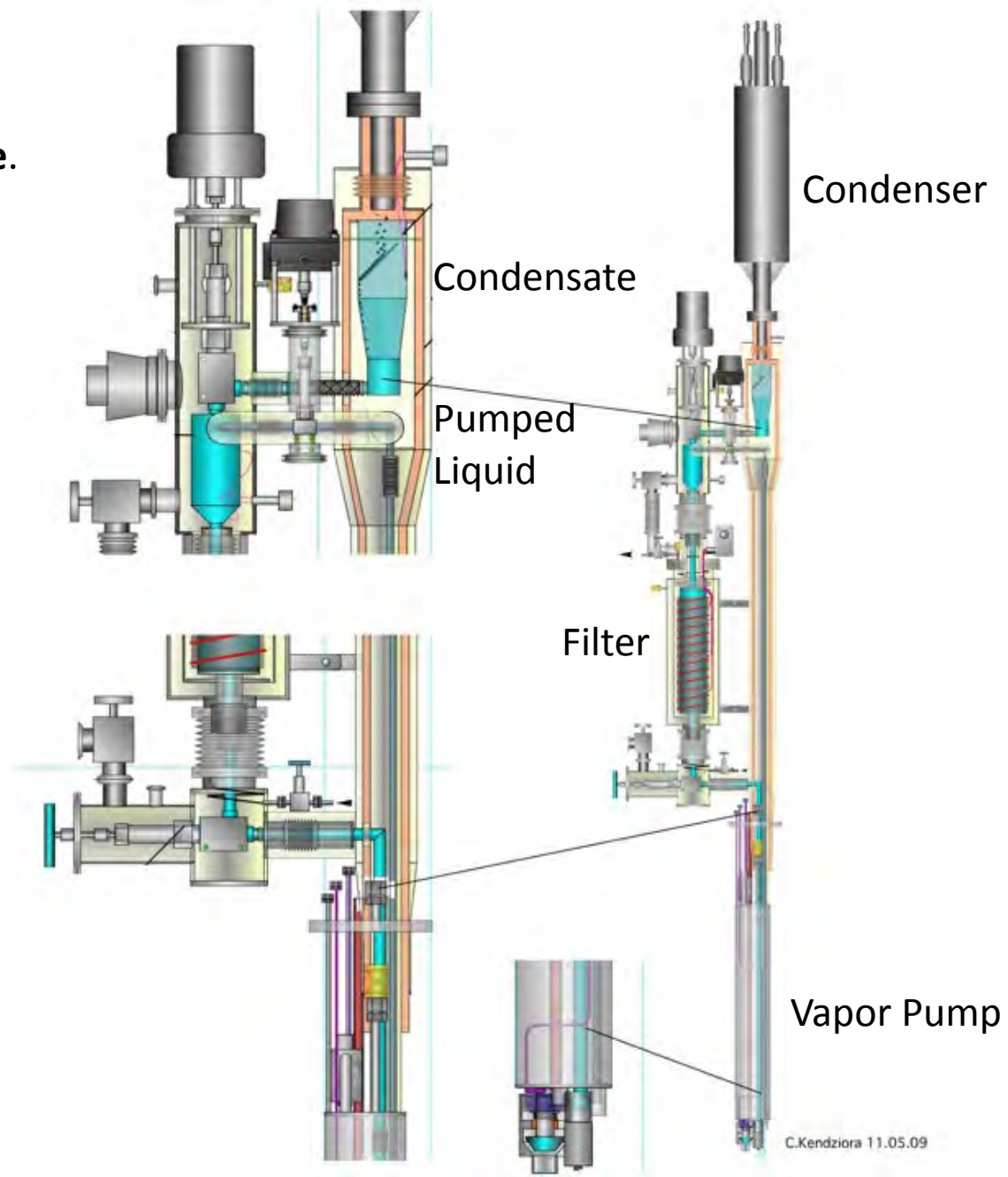


MicroBooNE

Single-wall, Foam Insulated Evacuatable Cryostat – 200 tons Liquid Argon

~ 9,000 Wire TPC – drift distance 2.6 m

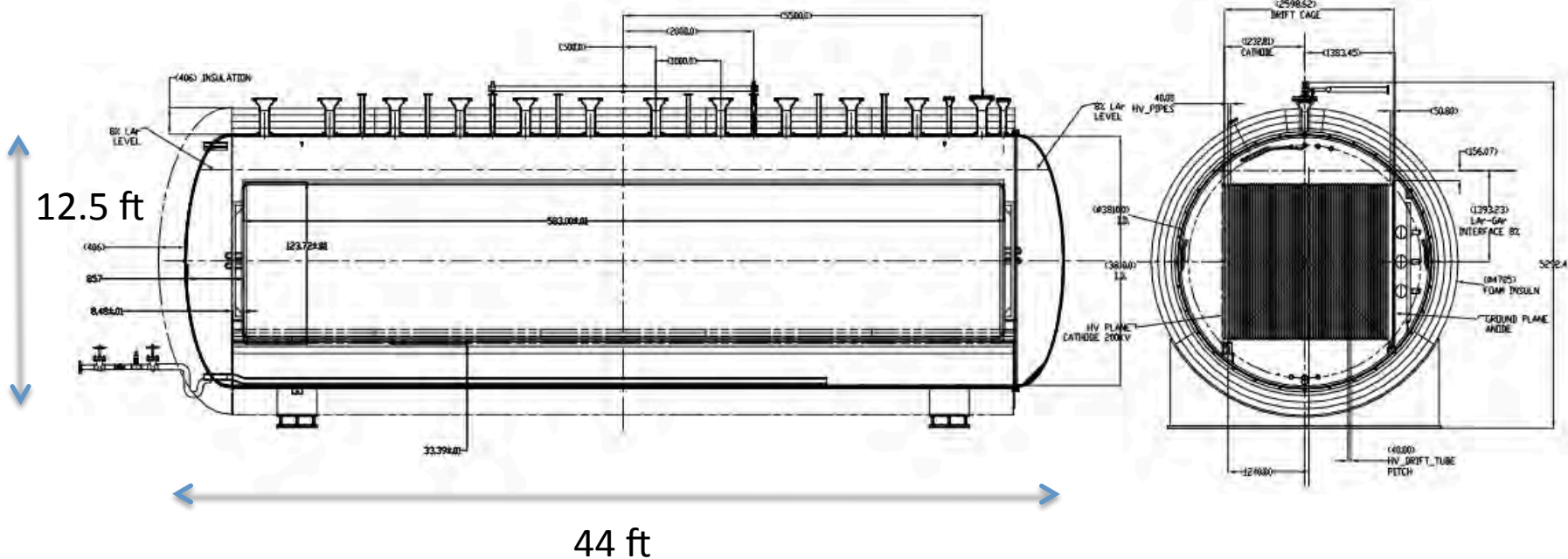

MicroBooNE Cryogenics and Purification Schematic

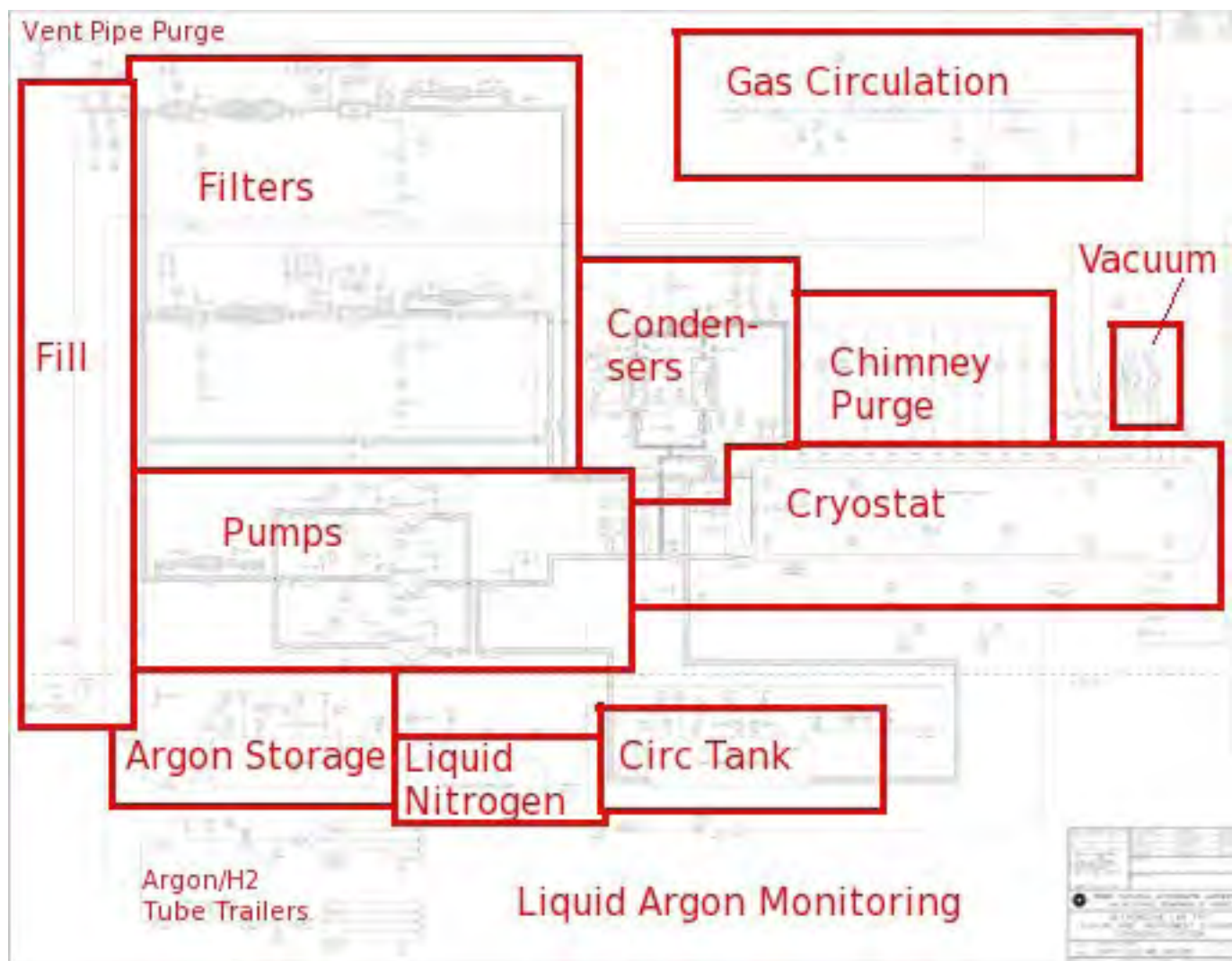

Role of MicroBooNE

Plan is to purge vessel to `purity' – details will be based on LAPD experience.

Liquid Recirculation and Purification specified - for 1 day cycle

Design of the Gas Circulation for Purification is under development – needs input from MTS and other purging studies.


If MicroBooNE purging to purity works, it is full validation of plan to go without evacuation

FIN

Purging Tests



Argon 'piston' well modeled – 4 ft rise /hr

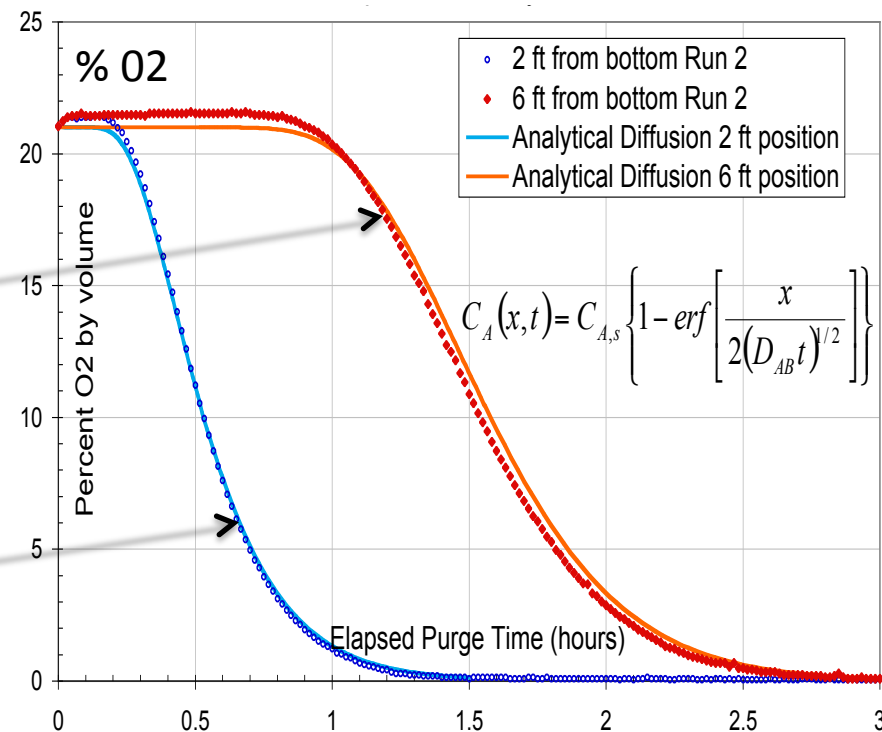

Outlet

Ambient Air



Inlet

O2 Monitors

Purge limit test of industrial tank

